

W E T - R I N G
V A C U U M S



CLASSIC

Introducing the Midmark ClassicSeries™ Wet-Ring Dental Vacuums

The Midmark ClassicSeries™ Wet-Ring Vacuums bring together the best features and most proven technologies into one uniform offering. These vacuums satisfy those who have grown to appreciate the high standards created by Apollo and Matrix brand wet-ring vacuums, offering:

- High Output, Proven Performance
- Durable Aluminum Bronze Construction
- Exclusive Internal Water Recycling
- Water Recycling By-Pass
- Stainless Steel Motor Shaft
- NFPA 99C Level 3 Compliant
- UL/CUL Classified

When it comes to selecting a wet-ring dental vacuum, Midmark is the natural choice. The ClassicSeries evolved from the Apollo and Matrix brand wet-ring vacuums, both of which have provided years of proven vacuum performance to thousands of satisfied dental professionals.

The ClassicSeries Wet-Ring Dental Vacuums come standard with many of the same value-driven features as their predecessors. Standard internal water recycling is an exclusive design that contributes to minimizing fresh water consumption – saving over three times the amount of fresh water consumed vs. non-water recycling vacuums. In addition, a valuable water recycler by-pass valve is utilized to help protect against the potentially harmful effects of water mineral deposits.


Aluminum bronze, rated harder than stainless steel (Brinell Hardness Scale), is used exclusively to ensure ClassicSeries wet-ring vacuums provide many years of durable operation. Stainless steel motor shafts and heavy-duty 1HP/2HP motors assure extra strength and reliability. Dual voltage (1HP vacuums) allows for versatility and ease of installation in all mechanical rooms.

Along with the best features brought forward from Apollo and Matrix, several new features elevate the ClassicSeries above all competitors. The component layout and design has been

simplified and re-organized to allow for extremely easy, convenient installation and service. From the strategically-designed left-hand/right-hand Twin vacuum installation layout to the convenient all-access front electrical panels, the ClassicSeries Wet-Ring Dental Vacuums definitely “raise the bar” in fit, finish, and serviceability.


You can depend on Midmark ClassicSeries Wet-Ring Dental Vacuums. Proven. Reliable. Trusted.

Operating Cost Over Time*

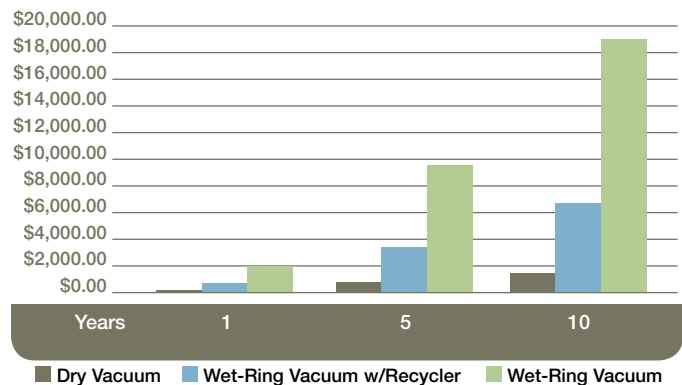

* 2,000 hours/yr; 57,600 gal/yr (Wet-Ring); 17,280 gal/yr (Wet-Ring w/Recycler)

Performance

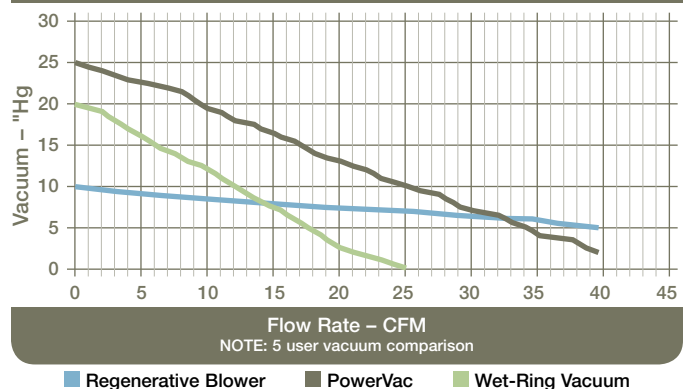



Redesigned intake manifold allows for smooth air flow and greater efficiency.



Separate Twin Model On/Off switches and Hour Meters provide operational versatility.



"Quick Connect" hose connections allow for tool-free serviceability. Water Recycler By-Pass Valve provides extra protection against water mineral deposits.

Specifications

Model	Max Users*	No. Vacuums	Total HP	Vacuum ("Hg)	Voltage	Amps (Per Vac)	Water Recycling	Overall Dimensions H x W x D	Shipping Weight (lbs.)
CV3	3	1	1	up to 12"	115/208-230	15/7.5	No	15" x 12" x 13"	70
CV3R	3	1	1	up to 12"	115/208-230	15/7.5	Yes	15" x 12" x 13"	70
CV5	5	1	2	up to 12"	208-230	7.5	No	17" x 12" x 13"	76
CV5R	5	1	2	up to 12"	208-230	7.5	Yes	17" x 12" x 13"	76
CV6	6	2	2	up to 12"	115/208-230	15/7.5	No	18.5" x 25.5" x 20"	188
CV6R	6	2	2	up to 12"	115/208-230	15/7.5	Yes	18.5" x 25.5" x 20"	188
CV10	10	2	4	up to 12"	208-230	7.5	No	18.5" x 25.5" x 20"	208
CV10R	10	2	4	up to 12"	208-230	7.5	Yes	18.5" x 25.5" x 20"	208

* HVE = 1 vacuum user; SE = .5 vacuum user; N₂O Scavenger = .5 vacuum user

Warranty

- 3 years/6,000-hour product warranty

Complies with NFPA 99C level 3



Midmark is an ISO 9001 Certified Company.

For more information or a demonstration, contact your Midmark dealer or Call: 1-800-MIDMARK Fax: 1-877-725-6495
Outside the U.S.A. Call: 1-937-526-3662 Fax: 1-937-526-8214
or visit our website at midmark.com

© 2007 Midmark Corporation
145 Mid County Drive Orchard Park, New York 14127-1726
Products subject to improvement changes without notice
Litho in U.S.A. 007-0652-00 Rev. A1 (4/07)



Because we care.