
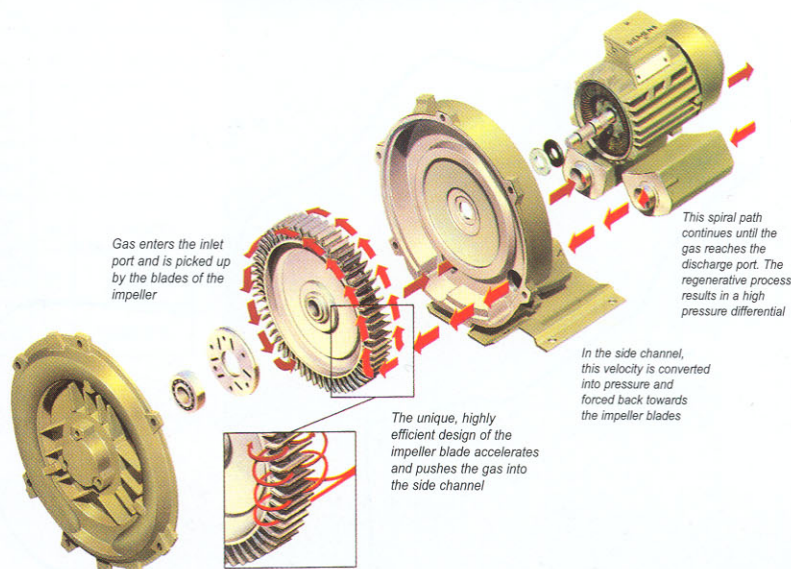
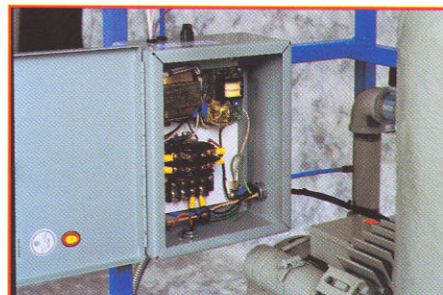
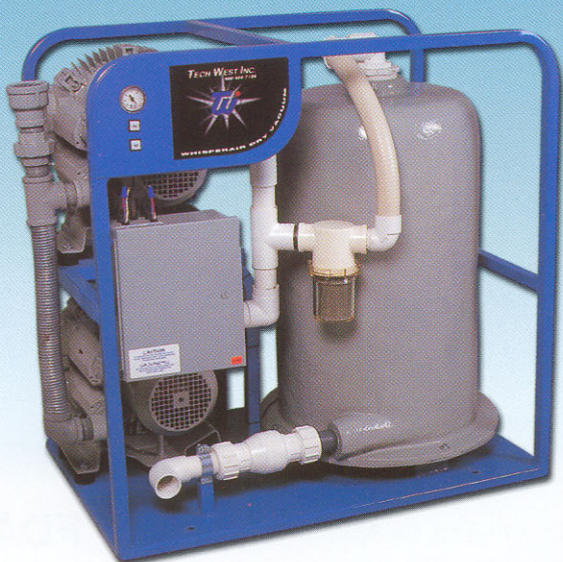


WHISPERAIR DRY VACUUM

- SINGLE 2-4 USER SYSTEMS
- DUAL 4-8 USER SYSTEMS
- QUIET BY DESIGN
- NO WATER USE
- NO LUBRICATION NEEDED
- NO VACUUM TRAP NEEDED

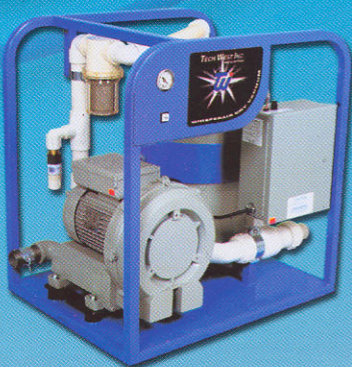


TW **TECH WEST INC.**

Manufacturers of Dental Vacuum
and Air Systems

2625 N. Argyle Ave. • Fresno, CA 93727
(559) 291-1650 • (800) 428-7139 • FAX (559) 348-9677

Tech West's Whisperair utilizes advanced technology in dry vacuum pump motors. Our system provides whisper quiet operation while supplying superior suction. The unit is completely self-contained and is shipped complete ready for installation. Our system can be installed on existing "wet" piping without the need to run new piping saving thousands of dollars. We use a solid state sensing system to drain the separator tank when it gets full. It is a simple and reliable design with no moving parts.

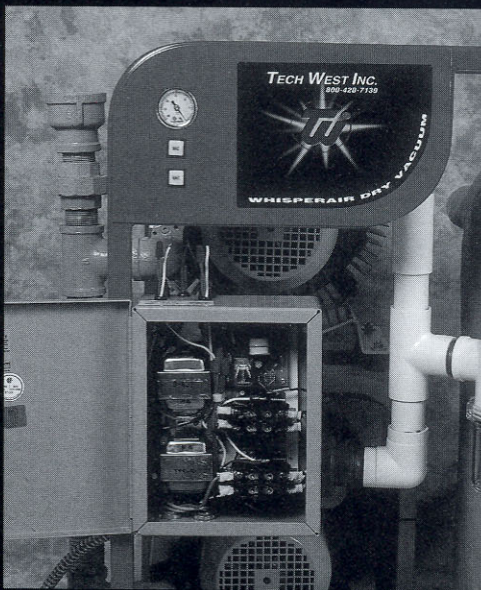

WHISPERAIR DRY VACUUM

FEATURES

Whisper quiet operation at 65 decibels.

The Whisperair is shipped completely assembled from the factory.

The Whisperair utilizes advanced dry vacuum motor technology without veins or lubrication.

The Whisperair doesn't use any water at all.

The Whisperair utilizes cyclonic separation of the evacuated waste and doesn't need a vacuum trap.

BENEFITS

This allows the equipment to be installed anywhere in the office without concern of noise.

This saves on installation time and costs.

This makes the Whisperair virtually maintenance free and will outlast other motor technologies.

This will save thousands of dollars in water and sewer costs over years.

This eliminates the need to constantly clean and replace the trap screens like other Dry-Vacs.

WHISPERAIR DRY VACUUM SYSTEMS

	MODEL NUMBER	# OF USERS	NO. OF MOTORS x H.P.	BREAKER SIZE	VOLTS	DB(A)	EQUIP. HEIGHT	EQUIP. DEPTH	EQUIP. WIDTH	EQUIP. WEIGHT	SHIPPING HEIGHT	SHIPPING DEPTH	SHIPPING WIDTH	SHIPPING WEIGHT
	VPD2SS	2	1 x 1HP	20	115/208/230	66	27	20	26	145	38	29	35	185
	VPD2S1	2	1 x 1HP	20	115	66	27	20	26	145	38	29	35	185
	VPD2S2	2	1 x 1HP	20	230	66	27	20	26	145	38	29	35	185
	VPD4D2	4	2 x 1HP	20	208/230	66	36	20	30	220	40	31	43	260
	VPD4S2	4	1 x 1.75HP	20	230	68	36	18	30	153	40	31	43	185
	VPD8D2	4	2 x 1.75HP	30	208/230	70	36	22	36	248	40	31	43	260

Revised 1/05

TECH WEST WISPERAIR DRY VACUUMS HAVE
A FULL 2 YEAR WARRANTY

Exceeds NFPA 99C Guidelines

TW **TECH WEST INC.**
Manufacturers of Dental Vacuum
and Air Systems

2625 N. Argyle Ave. • Fresno, CA 93727
(559) 291-1650 • (800) 428-7139 • FAX (559) 348-9677

FDA
RECOGNIZED